Remedy Massage System

Manual

VitaDora

LIVE. BETTER. LONGER.

Remedy Massage System

Congratulations on your purchase of the Vitadora Remedy Massage System.

The Vitadora Remedy Massage System provides an easy and effective way to stimulate your blood circulation, reduce pain and muscle soreness, decrease stress and improve your well-being. This compact and easy-to-use system provides the benefits of vibrational therapy combined with shiatsu massage and heat therapy, right in the comfort of your own home.


Caution: Please read all safety instructions carefully before use to ensure your safety.

- Keep this device, and all other electronic products, out of reach of children. Basic safety precautions must always be followed.
- Always unplug the device when not in use.
- Avoid overloading your home's electrical outlet.
- Keep both used and unused cords neat and secure to prevent damage.
- Do not use pins or other sharp objects on the massage system.
- Never use the Remedy Massage System if it is folded or crumpled.
- Do not stand on the unit.
- Keep electronic devices and power outlets away from water to prevent the risk of electric shock. Do not use the device when wet or damp. Do not use the Remedy Massager in humid or damp environments (e.g., in bathrooms or while showering).
- The Massage System is not intended to replace medical treatment.
- Consult your doctor before use if you are pregnant, have a pacemaker, have recently undergone surgery, suffer from blood clots, or have any other medical condition. If you experience any pain or discomfort during use, discontinue immediately and consult your doctor.
- Prolonged use may cause excessive heating, which could shorten the device's lifespan. After continuous use for 30 minutes, allow the device to cool for at least 30 minutes before using it again.
- If the temperature of the massager becomes too high for your comfort, stop use and allow the device to cool down.

These warnings align with the necessary safety guidelines to ensure safe usage of the product.

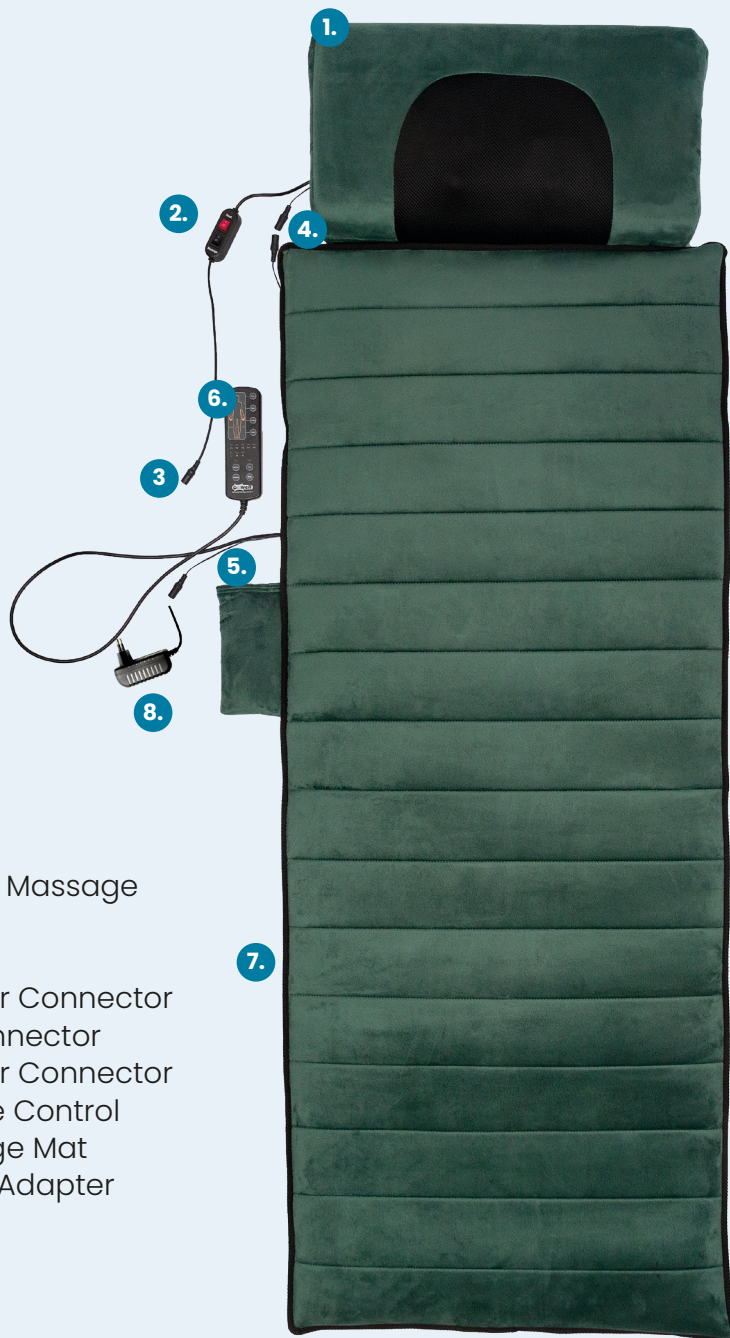1. Shiatsu Massage Pillow
2. Switch
3. Adapter Connector
4. Din Connector
5. Adapter Connector
6. Remote Control
7. Massage Mat
8. Socket Adapter


About the Remedy Massage System

The Vitadora Remedy Massage System contains two highly effective massage methods.

1. Massage mat with vibrational massage and heat therapy. The mat features strategically placed massage motors throughout the surface to specifically target every part of the body. It is possible to select a full body massage, or isolate certain muscle groups. There is a choice of 5 massage programs and 3 intensity levels.

2. Massage pillow with shiatsu massage and infrared heat. The 6 rolling massage heads provide instant relieve in your head and neck.

The Massage Pillow and massage mat can be used in combination with each other, but it is also possible to use them separately. The Massage Pillow can also be used for specific muscle groups, such as the back, legs and feet.

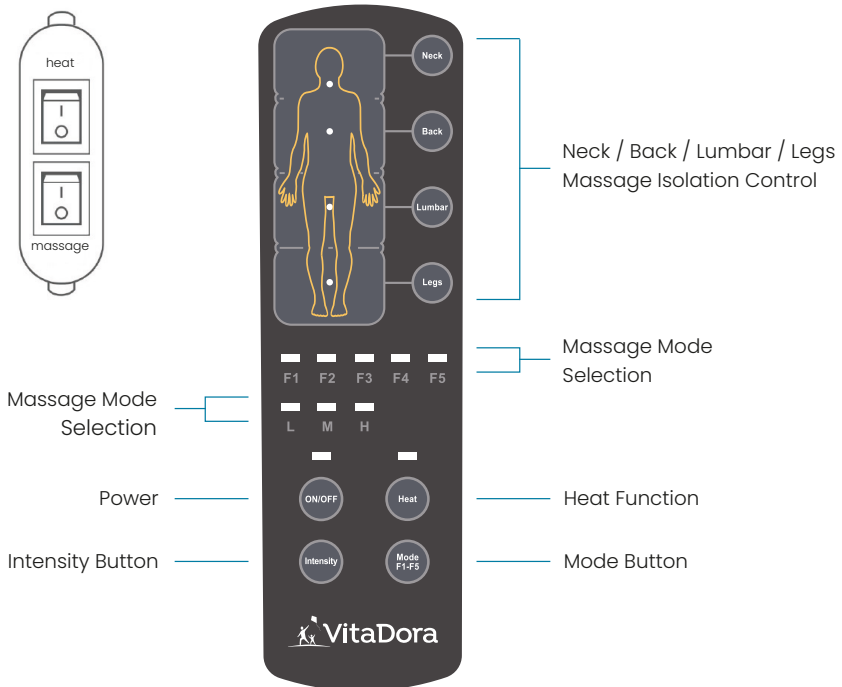Three Intensity Levels:

1. **L** for low intensity massages
2. **M** for medium intensity massages
3. **H** for high intensity massage

Five Massage Modes

1. Pulsate
2. Tapping
3. Rolling
4. Kneading
5. Intense

(After a few minutes the program will switch to Auto mode and automatically change through the different modes).

How to use

1. Unpack the Vitadora Remedy Massage System and all of its components.
2. Place the mat on a flat, stable surface such as the floor, a bed, sofa, or on a chair.
3. Plug the power adapter into the wall outlet and connect with the adapter connector of the Massage Pillow or the Massage Mat.
4. Turn on the device by pressing the ON/OFF button on the remote control.
5. Adjust the intensity, mode and optional heating on the remote control.
6. Lie down on the Remedy Massage System and position yourself comfortably. You can place a towel or a thin blanket between your body and the system to reduce the intensity of the vibration if necessary.
7. Relax and allow the vibrations to massage your body and neck. You can also choose to focus on specific areas such as your legs, lumbar, back and neck.
8. When you want to use the Massage Pillow or Massage Mat independently from each other, unzip the pillow from the mat, unplug the DIN Connector and connect the socket adapter with the device you like to use.
9. To turn off the device, press the ON/OFF button on the remote control.
10. When finished using the system, unplug the power cord and store the unit in a safe and dry place.

Manual Shiatsu - Nek- and shouldermassage


To fully enjoy the massage experience, it is important to pay attention to the position of your neck on the massage rollers.

Neck Position on the Massage Rollers

When using the Vitadora Remedy Massage System on your neck and shoulders, it is crucial to carefully choose the position of your neck on the massage rollers. Position your neck on the rollers so that it feels comfortable. If the shiatsu rollers are placed too high or too low on your neck, it may cause discomfort during the massage. We recommend trying different positions until you find a comfortable spot where the massage rollers reach the right areas of your neck and shoulders.

Adjusting the Position

If the position of the massage rollers doesn't feel perfect, do not hesitate to slightly reposition your body. Relaxation is key during the massage, so take your time to find the optimal position before starting.

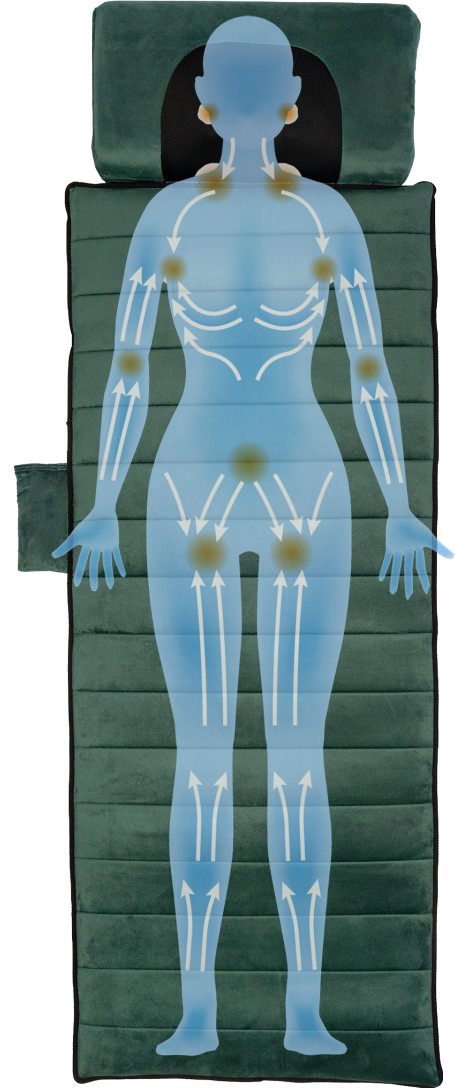Pressure Regulation

The Vitadora Remedy Massage System is designed to provide a deep shiatsu massage, but sometimes the pressure might feel too intense, especially if you are sensitive to firm massages. In such cases, we recommend placing a soft towel over the massage cushion.

Using a towel adds an extra layer of cushioning between the massage rollers and your neck and shoulders, reducing the pressure. You can adjust the thickness of the towel to suit your preference.

Personal Preferences

Everyone has their own preferences when it comes to massage intensity and comfort. The Vitadora Remedy Massage System is designed to adapt to your needs. Feel free to experiment with different positions and pressure levels to achieve the perfect massage experience for you.

Cleaning & Storing

Cleaning your unit

- To clean the outer shell of your massager, use a soft, slightly damp sponge.
- Do not allow water or any liquid to come in contact with the unit.
- DO NOT use any chemicals or harsh abrasives such as cleaners and brushes to clean.
- NEVER immerse your Remedy Massage System or any other electronic device in water.
- Never clean when the unit is turned on.
- Make sure the massager is turned off and unplugged from the adapter before cleaning.

Storing you unit

- Unplug the cord from the power source.
- Loosely fold the massage mat.
- Place the massager, along with the cord, in its carry bag, or in a safe, dry and cool place.
- AVOID contact with sharp edges or pointed objects.

DISPOSAL

Washer, accessories and packaging should be sorted for environmentally friendly recycling. This product contains a lithium-ion battery. Do not dispose of battery packs into household waste, water or fire. Please dispose of them at a public recycling centre.

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