motionsolution

full spine stretch & mobilization mat

Dr. Michael Ho

Spine Care

The Solution for a
Fit & Flexible Spine

User manual | P. 2

Bruksanvisning | P. 13

Brukermanual | P. 24

Brugervejledning | P. 34

Käyttöopas | P. 46

Dear Friends,

Thank you for choosing my DR-HO'S MotionSolution Mat. I truly appreciate your trust in me.

I developed this stretch and relaxation mat based on over 30 years of clinical experience treating patients with sore and stiff spines. You'll love how easy it is to use, and you'll feel incredible comfort and relief both during and after using the mat.

DR-HO'S MotionSolution Mat is designed to relieve tension and discomfort from your head, neck, upper shoulders, mid back, lower back, and hip areas. It offers one of the most comfortable and enjoyable solutions for a stiff and sore spine.

The gentle mobilizing motion of DR-HO'S mat helps stretch and relax your spinal muscles, joints, and soft tissues, relieving tension, stiffness, and soreness exactly where you need it.

For your head, neck, and upper

shoulder areas: The mat's flexion and extension, side-to-side rotation, and progressive bouncing and shaking motions will work wonders to release tension, soreness, and stiffness in your neck and shoulders with every movement.

For your upper back and posture:

The mat will twist, lift, and extend your upper and middle back, reducing slouching and addressing any hump that may have formed in your upper back. After the treatment, you'll feel great and find it easier to stand tall.

For your lower back and hip area:

With the heated back support roll positioned under the arch of your back, your muscles will warm up and loosen. The mat will then gradually lift your back with gentle pauses and bounces, providing an incredible stretch and extension. It will also push and pull on your hips, twisting and rotating your lower back

to loosen the muscles, joints, and soft tissues that cause stiffness and soreness in your back and hips. The combination of heat therapy, progressive traction, and passive stretching will feel amazing, especially for anyone dealing with chronic back or disc issues.

I encourage you to use DR-HO'S MotionSolution Mat regularly to maintain a healthy range of motion and promote spinal health. For best comfort, place the mat on a flat surface such as a bed, couch, or guestroom mattress. You can also use it on the floor if that feels comfortable. I suggest using it before bedtime and again in the morning to keep your spine feeling great throughout the day. Each session can last 20 to 30 minutes, or even longer if you like. It's drug-free and can be used daily without any side effects. You'll feel so relaxed while using the mat, and after each treatment, you'll stand tall with good posture, feel loose and agile, and be ready to move with ease.

With DR-HO'S MotionSolution Mat, you can stretch and relax your spine anytime and enjoy life to the fullest!

Warm regards,
Dr. Michael Ho

P.S. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact DR-HO'S support team with your product model number.

IMPORTANT SAFEGUARDS

When using your DR-HO'S MotionSolution® Mat, always observe basic safety precautions.

Warnings

- Avoid exposing the product to water, other liquids, direct sunlight, or high temperatures.
 - Never leave MotionSolution® unattended while it is in use.
 - Use the product only for its intended purposes.
 - This system shall not be used or played with by children. Keep the system and its cord out of reach of children.
1. Use of the system in locations with different electrical specifications may result in damage to the product so please check specifications before operation. Do not use this product with any type of adapter or voltage converter device; this is considered unauthorized modification of the product and as such voids the warranty.
 2. Do not pull, twist, or damage the power cord and don't allow to touch hot surfaces. Periodically inspect cord and plug. Do not operate with a damaged cord or plug. A damaged cord or plug must be replaced by a qualified person in order to avoid a hazard.
 3. If the system malfunctions, or is damaged in any manner, contact

Customer Service for assistance.

Have the product repaired only by a qualified professional if it becomes defective. Regularly inspect the product for signs of wear or damage before use.


4. Do not use the mat outdoors or in inclement weather, and always operate on a level surface.
5. This system may be used, cleaned and/or maintained by persons with reduced physical, sensory or mental abilities or lack of experience and/or knowledge, if they have been given supervision or instruction concerning use of the system in a safe way and if they understand the hazards involved.
6. Never stand on the air mat. Lie down on the mat as described later in this manual.
7. Switch OFF the MotionSolution and remove the plug from the socket every time you finish using the device and before cleaning it. Proceed as follows: take hold of the plug (not the cable) and pull the plug from the socket. Never pull on the power cable.
8. Contact your supplier in case of problems.

CONTENTS

- 1 x DR-HO'S MotionSolution Air Massage mat
- 1 x Adapter
- 1 x User manual

PRODUCT DIAGRAM

1. Neck Mobilizer
2. Heated Back Support
3. Back Mobilizer
4. Blue Buttock zones
5. Leg rest zones
6. Remote Control
7. Power cord with adapter
8. Zipper, for easy storage

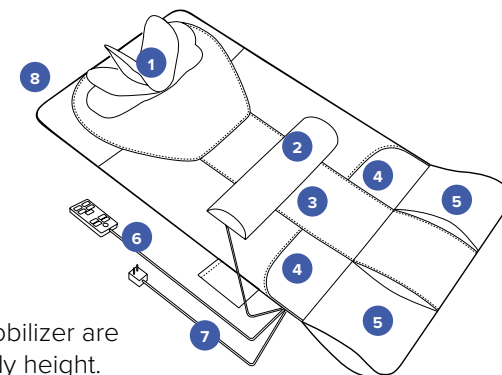The Neck Mobilizer and Heated Back Mobilizer are both adjustable in height to suit your body height. The Heated Back Mobilizer can also be removed entirely, for example if you prefer lying flat or want to store it separately.

The raised areas on the mat provide thigh support but do not create any motion; they house the pumps and offer a comfortable resting place for your legs.

PRODUCT SPECIFICATIONS

Product code Model	H08646 MSM-8801
Name	DR-HO'S MotionSolution Mat
Voltage	12V  2.5A 30W
Default timer	10±1 minute
Functions	Specific treatment for Neck / Shoulders, Mid Back and-Lower back, plus 4 selectable Auto Modulated Motion Programs
Materials	Vegan Leather, Suede, Wool
Dimensions	120 x 60 cm
Weight	4,68 kg

PRECAUTIONS AND CONTRAINDICATIONS

- Do not use the MotionSolution Mat if you have fractured bones, are post-surgery, have skin ulcerations or other health conditions that are not suitable for heat or spinal movement exercise.
- If you feel any discomfort during or after using the system, stop and rest for a day or two before using it again. If the discomfort continues, stop using the system.
- As always, please consult your doctor regarding any health issues or concerns before starting new exercises, using this mat, or trying any new devices.

Disclaimer:

DR-HO'S products are not intended to diagnose, treat or cure diseases and disorders. In the event of health problems, consult your doctor for proper diagnosis and treatment.

WHAT TO EXPECT WHEN YOU START

If you have significant stiffness in your neck or lower back, using the mat to stretch these areas may cause some initial discomfort. Don't be concerned, this is a normal reaction when starting any new form of exercise. The discomfort is a result of renewed movement in stiff joints and tight muscles, and it's actually a good sign that your body is beginning to respond.

With regular use of the mat for 1 to 2 weeks, you should start to feel looser and more comfortable, without the initial discomfort. However, if the discomfort persists or if you have any health concerns, it's best to consult your doctor before continuing use.

Tips for a good start

To help your body gently get used to the mat:

- Use the mat for 10 minutes per session, once or twice a day, during the first 2 weeks.
- Position the Neck Mobilizer under your neck, not under your head.
- Use the mat without the red Back Mobilizer during the first 2 weeks.

Note: You can easily remove the red Back Mobilizer by loosening the Velcro strap at the bottom. Don't forget to unplug the wire connected to the back support.

HOW TO POSITION THE MAT AND YOURSELF

For your comfort, place the mat on a flat surface where you can lie down easily, a bed, couch, or guestroom mattress works well. Some people also prefer placing it on the floor, especially if they're used to yoga or floor exercises. Ideally, choose a spot where the mat can stay in place for regular use, but above all, choose whatever is most convenient for you at home.

Cable connections

- The remote control (RC) is already connected.
- Connect the matching colored connectors together.

Next steps:

1. Position your buttocks over the two blue zones.
2. Lie down with your knees bent or your legs straight out.
3. Ensure the Back Support (the red adjustable bar with heat function) is positioned under the arch of your lower back by sliding it up or down.

4. Ensure the Neck Mobilizer supports your neck by sliding it up or down to suit your height.


Note: Experiment with different positions to find the most comfortable alignment for your spine.

The best way to get off the mat (and out of your bed):

1. Turn onto your side, allowing your legs to hang over the edge of the bed or couch.
2. Push up using your elbow and your other hand.

This method is easy and is the best for your neck and lower back!

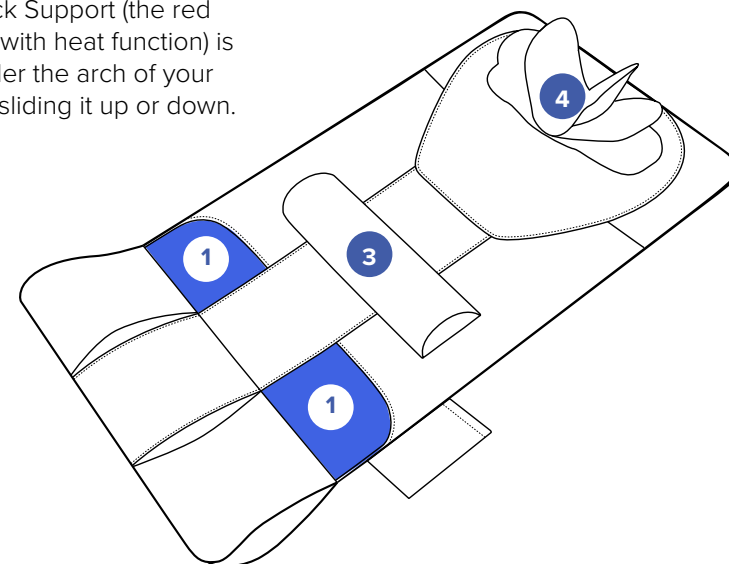INTENDED USE

DR-HO'S MotionSolution Mat is a passive motion stretch and relaxation mat designed to stretch and relax the muscles, joints, and soft tissues along your spine. It works to loosen your spine, helping to keep it feeling good and comfortable.

The mat is engineered with a series of Body Mobilizers, precisely positioned to automatically, sequentially, and continuously move specific areas of your spine in six planes of motion. These include flexion and extension, right and left rotation, and right and left twisting, along with something you've likely never experienced before: Progressive Elevation Extension with a gentle bouncing motion. It feels amazing to have your spine elongated, passively moved, and stretched in all these different directions.

Invented by Dr. Michael Ho, a Doctor of Chiropractic Medicine, this mat is programmed to move and stretch the spine in a way similar to how a professional Spine Therapist would.

Additionally, DR-HO'S MotionSolution Mat includes a special heated back support zone that warms and relaxes the muscles in the lower back. Simultaneously, it lifts the lower back with a unique Progressive Traction Elevation motion.

The combination of heat and traction works to:

- stretch tight muscles,
- open stiff joints,
- restore elasticity in ligaments and tendons, and
- decompress spinal discs.

This makes your back feel relaxed and comfortable, contributing to better physical health, enhanced athletic performance, and improved quality of life by enabling easier movement and reducing discomfort.

Advice:

It is recommended to use the MotionSolution Mat 20 to 30 minutes each session, and longer if desired.

STEP BY STEP INSTRUCTIONS

You can use the mat with the default settings, or you can select specific segments of the spine to treat. Alternatively, you can choose one of the four auto-modulated programs: P-1, P-2, P-3, or P-4.

There are three levels of intensity and options for treatment durations of 10, 20, or 30 minutes, all with an auto-shutoff feature.

First, press this button to turn the power on. The mat will enter standby mode with the default settings set to Program P1, medium intensity, and a treatment time of 10 minutes.

To start treatment with these default settings, press this button a second time.

Alternatively, press the buttons below for customized treatment:

Press this button if you want to select a specific area of the spine to be treated. The corresponding area will light up for one of the following regions:

- Neck & Shoulders,
- Upper & Mid Back,
- Lower Back & Hips.

Press this button to select the treatment time of 10, 20, or 30 minutes.

Press one of the 4 buttons to select a specific program: P1, P2, P3 or P4:

Body Extension Motion

Benefit from the auto-modulated mode 'Body Extension Motion', which focuses on extending the spine to help relieve tension, stiffness, and soreness.

P2 Body Waving Motion

The auto-modulated mode 'Body Waving Motion' focuses on restoring the joint movements of the spine, relieving the joints and muscles, and helping to restore flexibility and agility.

P3 Body Twisting Motion

Experience the benefits of the auto-modulated mode 'Body Twisting Motion,' which focuses on stretching and releasing tight ligaments around the spine, helping to restore elasticity and mobility.

P4 Body Energizing Motion

Discover the benefits of the auto-modulated mode 'Body Energizing Motion,' which focuses on bouncing and shaking movements to help loosen tight fibrous tissues from old injuries and aging. This mode also supports good postural alignment and balances all range of movements along the spine.

 Press this button to pause the operation of the mat. Press it again to resume.

 Press this button to select the intensity level: L1, L2, or L3. Choose the level that makes your treatment as comfortable as you wish, allowing you to customize your experience according to your needs.

 Press this button to activate the soothing heat for the back support. The warmth not only enhances your comfort but also makes your treatment even more effective. If you feel you've had enough, just press it again to turn off the heat.

PROBLEM SOLVING

Problem	Cause	Solution
Abnormal smell.	New foam and cover materials may have a slight odor when first opened, but this is completely normal.	To help reduce the odor, air out the mat in an open room for 24 hours, and it should soon be freshened up.
Obvious creases on the mat.	These will be caused by packing and rolling it up for shipping.	Simply unfold the mat and lay it flat to smooth it out.
The surface of the mesh fabric appears uneven after use.	This could be due to trapped air in the Body Mobilizers.	To fix this, gently press down on the mat to help release any air pockets.
When I started using the mat, the movements felt unfamiliar and a bit strange.	Your spine is being gently mobilized in ways that support flexibility and relieve tension. This can feel unusual at first, but it's a sign your body is being guided through healthy motion.	Begin with shorter sessions to let your body ease in. As your spine gets used to the movements, they'll feel more natural and deeply relaxing. You may notice improved posture, more flexibility, and a refreshed, upright feeling after each use.
The mat stops working automatically.	Please note that the mat is designed to shut off automatically. It will turn off after 10, 20, or 30 minutes of operation, depending on the timer you've selected.	If you'd like to continue your session, just press the ON/OFF button to restart the mat.

MAINTENANCE AND CLEANING

- Store the mat in a safe, dry, and well-ventilated space.
- Use the zipper to fold the mat into three sections and zip up the sides for easy storage and travel.
- If you won't be using the mat for an extended period, cover it with a dust-proof cloth or zip it up and store it in its original packaging to prevent dust accumulation.
- Always disconnect the mat from the power supply before cleaning.
- Avoid using corrosive detergents; instead, use a damp cloth for cleaning.

WARRANTY AND COMMUNICATIONS

Warranty: one year against factory and material defects. The warranty does not cover faults caused by incorrect handling, damage, assembly or attempted repairs. The same applies to normal wear and tear. This warranty does not affect applicable legal rights.

For questions, assistance, support, or warranty claims, please contact us via the email below.

VitaDora UK Ltd.
www.drho.co.uk
E: support@drho.co.uk

Happy with your DR-HO'S MotionSolution Mat®? Share your posts, tips & tricks! Like & follow us @drho.uk & tag your posts

#MyDrHo #DrHoLovers #DrHoNow

#myflexiblespine #MotionSolution #SpineCare

SAVE THESE INSTRUCTIONS FOR HOUSEHOLD USE ONLY

THIS MARKING INDICATES THAT THIS PRODUCT SHOULD NOT BE DISPOSED WITH OTHER HOUSEHOLD WASTES THROUGHOUT THE EU. TO PREVENT POSSIBLE HARM TO THE ENVIRONMENT OR HUMAN HEALTH FROM UNCONTROLLED WASTE DISPOSAL, RECYCLE IT RESPONSIBLY TO PROMOTE THE SUSTAINABLE REUSE OF MATERIAL RESOURCES. TO RETURN YOUR USED DEVICE. PLEASE USE THE RETURN AND COLLECTION SYSTEMS OR CONTACT THE RETAILER WHERE THE PRODUCT WAS PURCHASED. THEY WILL TAKE CARE FOR ENVIRONMENTAL SAFE RECYCLING.

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VitaDora (Myori Nordic AB)
Suppl.: DRTV Asia Limited, 9/F, Flat/Room